
Occlutech ASD Occluder

Product specifications





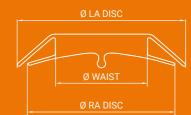

Article No.	Defect Size	Ø Waist	Ø LA Disc	Ø RA Disc	Sheath Size*	Pusher Article No. Pistol Pusher	(color)** Flex Pusher II
INO.						Pistoi Pushei	riex Pustiei II
29ASD04	D≤4	4	11		7	55PP090 🔵	51FP060 🔵
29ASD05	4 <d≤5< td=""><td>5</td><td>14</td><td>11</td><td>7</td><td>55PP090 🔵</td><td>51FP060 🔵</td></d≤5<>	5	14	11	7	55PP090 🔵	51FP060 🔵
29ASD06	5 <d≤6< td=""><td></td><td>16.5</td><td>12.5</td><td>7</td><td>55PP125 🔵</td><td>51FP100 🔵</td></d≤6<>		16.5	12.5	7	55PP125 🔵	51FP100 🔵
29ASD07	6 <d≤7.5< td=""><td>7.5</td><td>18</td><td>14</td><td>7</td><td>55PP125 🔵</td><td>51FP100 🔵</td></d≤7.5<>	7.5	18	14	7	55PP125 🔵	51FP100 🔵
29ASD09	7.5 <d≤9< td=""><td></td><td>20.5</td><td>16.5</td><td>7</td><td>55PP125 🔵</td><td>51FP100 🔵</td></d≤9<>		20.5	16.5	7	55PP125 🔵	51FP100 🔵
29ASD10	9 <d≤10.5< td=""><td>10.5</td><td>22</td><td>18</td><td>7</td><td>55PP125 🔵</td><td>51FP100 🔵</td></d≤10.5<>	10.5	22	18	7	55PP125 🔵	51FP100 🔵
29ASD12	10.5 <d≤12< td=""><td>12</td><td>27</td><td>23</td><td></td><td>55PP165 🔾</td><td>51FP120 🔘</td></d≤12<>	12	27	23		55PP165 🔾	51FP120 🔘
29ASD13	12 <d≤13.5< td=""><td>13.5</td><td>28.5</td><td>24.5</td><td>9</td><td>55PP165 🔘</td><td>51FP120 🔘</td></d≤13.5<>	13.5	28.5	24.5	9	55PP165 🔘	51FP120 🔘
29ASD15	12 <d≤15< td=""><td>15</td><td>30</td><td>26</td><td></td><td>55PP165 O</td><td>51FP120 🔘</td></d≤15<>	15	30	26		55PP165 O	51FP120 🔘
29ASD16	15 <d≤16.5< td=""><td>16.5</td><td>31.5</td><td>27.5</td><td>9</td><td>55PP165 O</td><td>51FP120 🔘</td></d≤16.5<>	16.5	31.5	27.5	9	55PP165 O	51FP120 🔘
29ASD18	15 <d≤18< td=""><td>18</td><td>33</td><td>29</td><td></td><td>55PP165 O</td><td>51FP120 O</td></d≤18<>	18	33	29		55PP165 O	51FP120 O
29ASD19	16.5 <d≤19.5< td=""><td>19.5</td><td>34.5</td><td>30.5</td><td>10</td><td>55PP165 🔘</td><td>51FP120 🔘</td></d≤19.5<>	19.5	34.5	30.5	10	55PP165 🔘	51FP120 🔘
29ASD21	18 <d≤21< td=""><td>21</td><td>36</td><td>32</td><td>11</td><td>55PP185 🔘</td><td>51FP150 🔘</td></d≤21<>	21	36	32	11	55PP185 🔘	51FP150 🔘
29ASD24	21 <d≤24< td=""><td>24</td><td>39</td><td>35</td><td>11</td><td>55PP185 O</td><td>51FP150 🔘</td></d≤24<>	24	39	35	11	55PP185 O	51FP150 🔘
29ASD27	24 <d≤27< td=""><td>27</td><td>42</td><td>38</td><td>12</td><td>55PP210 🔵</td><td>51FP160 🔵</td></d≤27<>	27	42	38	12	55PP210 🔵	51FP160 🔵
29ASD30	27 <d≤30< td=""><td>30</td><td>45</td><td>41</td><td>12</td><td>55PP210 🔵</td><td>51FP160 🔘</td></d≤30<>	30	45	41	12	55PP210 🔵	51FP160 🔘
29ASD33	30 <d≤33< td=""><td>33</td><td>48</td><td>43</td><td>12</td><td>55PP210 🔵</td><td>51FP160 🔵</td></d≤33<>	33	48	43	12	55PP210 🔵	51FP160 🔵
29ASD36	33 <d≤36< td=""><td>36</td><td>52</td><td>46</td><td>12</td><td>55PP210 🔵</td><td>51FP160 🔵</td></d≤36<>	36	52	46	12	55PP210 🔵	51FP160 🔵
29ASD39	36 <d≤39< td=""><td>39</td><td>54</td><td>49</td><td>12</td><td>55PP210 🔵</td><td>51FP160 🔵</td></d≤39<>	39	54	49	12	55PP210 🔵	51FP160 🔵
29ASD40	39 <d≤40< td=""><td>40</td><td>55</td><td>50</td><td>12</td><td>55PP210 🔵</td><td>51FP160 🔵</td></d≤40<>	40	55	50	12	55PP210 🔵	51FP160 🔵

Quality made in Europe

Occlutech adheres to the highest quality standards. We never compromise on quality. Our dedication to patients and the medical community is the main driver of our vision. We want to be at the forefront of developing technologies and products that matter, that make a difference and that help our customers to be at the leading edge of better treatment for structural heart disease.





Perfecting Performance



Occlutech ASD Occluder




Designed for safety and excellence



The design of the Flex II ASD Occluder results in an ideal septum alignment which increases its feasibility and patient safety during implantation.¹

- · Optimized braiding pattern and less material at left atrial disc.
- Titanium oxide-covered Nitinol shows the best values of biocompatibility² and minimum Nickel release.³
- Spunbonded PET-patch allows for fast endothelialization.
- The unique ball-connection between pusher and Occluder safely locks it, while it freely follows the anatomy. Once in place it is easy and fast to deliver.

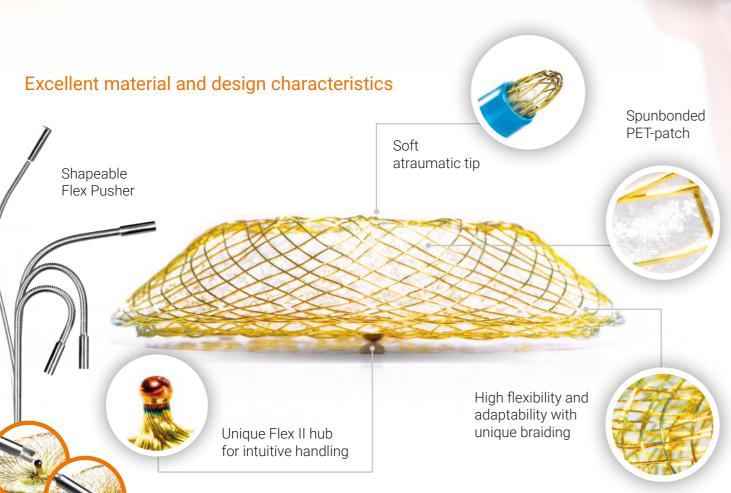







Recent studies found that transcatheter closure of isolated secundum atrial septal defects using the novel design of the Occlutech ASD Occluder was safe, effective and had an excellent outcome during the 6 months follow-up period.^{4,5}


- · Lower pull in forces and improved delivery sheath compatibility.
- · Reduced material exposure.
- Significantly reduced mean procedure and fluoroscopic time compared to other devices (34%).⁶

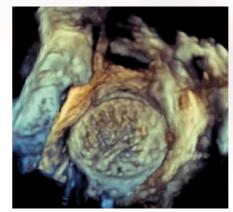

3D TEE image of an implanted Occlutech ASD Occluder.

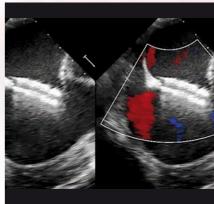

Perfect alignment of the Occluder to the

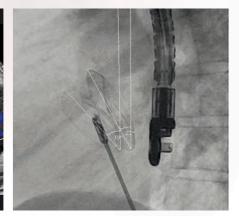

50° angulation of Flex Pusher vs Occluder hub.

